

12/13/02

USER MANUAL

MP-3R

GEMINI



merits

WELCOME

Welcome aboard the MP3R. We believe that you have chosen a product that will enrich your life by allowing you the freedom to go places you have always wished to go. This unique user friendly MP3R is designed for maximum maneuverability, even in narrow spaces. With your new MP3R, you will enjoy the freedom of traveling and mobility.

 **merits**

INDEX

I.	SAFETY INSTRUCTION	...	1
II.	EMI	...	4
III.	ENVIRONMENTAL CONDITIONS	...	5
IV.	SPECIFICATION	...	6
V.	ASSEMBLY & DISASSEMBLY		
	A). Install Battery & Ignition key Connection	...	7
	B). Install Seat	...	8
	C). Insert the Height and Width adjustable Armrests	...	8
	D). Install the Controller	...	9
VI.	ADJUSTMENTS FOR SEATING CONFORT		
	A). Seat Height	...	11
	B). Armrest Width, Height, and Angle Adjustment	...	11
	C). Foot-Platform	...	13
	D). Controller	...	14
	E). Backrest	...	14
	F). Anti-tipper	...	14
VII.	OPERATING		
	A). Driving	...	16
	B). Controller Display	...	17
	C). Freewheeling	...	19
	D). Thermal Rollback	...	19
	E). Main Circuit Breaker	...	20
	F). Seat Rotation	...	20

VIII. BATTERY and CHARGER

A). General Battery Information	... 21
B). Character of this type Battery	... 21
C). Portable Charging Instructions	... 23

IX. ACCESSORIES

... 24

X. MAINTENANCE and REPAIR

A). Maintenance	... 25
B). Repair	... 25

SAFETY INSTRUCTIONS

Your **Merits** MP3R is a battery-operated personal mobility vehicle. Therefore, it is required to exercise caution and consideration when operating it as to ensure your own personal safety and that of others around you.



The following are the rules for safe operation of your **Merits** Gemini:

1. Make sure that the unit is off before getting in or out of it. This will eliminate the possibility of accidentally activating the joystick and causing injury to yourself or others.
2. Always use a seat-belt and keep your arms on or inside the armrests and feet on the footrest at all times.
3. Set the speed control knob according to your driving ability and the environment where you are going to operate it in. We recommend that you keep the speed at the slowest position (fully counter-clockwise) until you are familiar with the driving characteristics of this vehicle. We also recommend that you do so when using your Merits MP3R indoors.
4. Always climb or descend gradients perpendicular to the slope. Never drive across the slope.
5. Do not attempt to climb or descend gradients greater than 12 (slopes of 1/5).
6. Do not carry any passengers under any circumstance. Otherwise, your MP3R will become unstable. Personal injury may result.
7. Never drive on the roadway.
8. Always cross street intersections via the most direct route making sure that you are visible to motor traffic.
9. Leave and join curb-cuts perpendicular to the road.
10. Use caution when driving over soft, uneven or unprotected surfaces such as grass, gravel and decks.

11. Never occupy your **MERITS MP3R** when transporting it in a motor vehicle. Make sure that it is securely strapped with an approved tie-down system.
12. Do not operate your vehicle if it isn't functioning properly.
13. Do not connect, or allow someone else to connect, any electrical or mechanical device to the vehicle. Failure to do this may result in injury and will void the warranty.
14. Never use electronic radio transmitters such as CB's, walkie-talkies, portable computers or cellular phones while using the vehicle without first turning the vehicle off.
15. Check with your physician if you are taking any medication that may affect your adjustment and ability to operate your Gemini.

Electromagnetic Interference (EMI) from Radio Wave Sources

The rapid development of electronics, especially in the area of communications, has saturated our environment with electromagnetic (radio) waves that are emitted by television, radio and communication signals. These EM waves are invisible and their strength increases as one approaches the source. All electrical conductors act as antennas to the EM signals and, to varying degrees, all power wheelchairs and scooters are susceptible to electromagnetic interference (EMI). This interference could result in abnormal, unintentional movement and/or erratic control of the vehicle. The United States Food and Drug Administration (FDA) suggests that the following statement be incorporated to the user's manual for all power wheelchairs like the Merits GEMINI:

Powered wheelchairs and motorized scooters (in this text, both will be referred to as powered wheelchairs) may be susceptible to electromagnetic interference (EMI), which is interfering electromagnetic energy emitted from sources such as radio stations, TV stations, amateur radio (HAM) transmitters, two-way radios and cellular phones. The interference (from radio wave sources) can cause the powered wheelchair to release its brakes, move by itself or move in unintended directions. It can also permanently damage the powered wheelchair's

control system. The intensity of the EM energy can be measured in volts per meter (V/m). Each powered wheelchair can resist EMI up to a certain intensity. This is called the "immunity level." The higher the immunity level, the greater the protection. At this time, current technology is capable of providing at least 20 V/m of immunity level, which would provide useful protection against common sources of radiated EMI.

Following the warnings listed below should reduce the chance of unintended brake release or powered wheelchair movement that could result in serious injury:

- 1) Do not turn on hand-held personal communication devices such as citizens band (CB) radios and cellular phones while the powered wheelchair is turned on.
- 2) Be aware of nearby transmitters such as radio or TV stations and try to avoid coming close to them.
- 3) If unintended movement or brake release occurs, turn the powered wheelchair off as soon as it is safe.
- 4) Be aware that adding accessories or components, or modifying the powered wheelchair, may make it more susceptible to interference from radio wave sources. (Note: there is no easy way to evaluate their effect on the overall immunity of the powered wheelchair).
- 5) Report all incidents of unintended movement or brake release to the powered wheelchair manufacturer, and note whether there is a radio wave source nearby.

TURN OFF YOUR GEMINI AS SOON AS POSSIBLE WHEN EXPERIENCING FOLLOWS:

1. Unintentional motions.
2. Unintended or uncontrollable direction.
3. Unexpected brake release.

FDA has written to the manufactures of power wheelchairs, asking them to test their new products to be sure they provide a reasonable degree of immunity against EMI. The letter says that powered wheelchair should have an immunity level of at least 20 V/m, which provide a reasonable degree of protection against the more common sources of EMI. The higher the level, the greater the protection.

Your MP3 R has an immunity level of 20 V/m which should be a very good of protection against the EMI.

ENVIRONMENTAL CONDITIONS

The environmental conditions may affect the safety and performance of your **MP3R**. Water and extreme temperatures are the main elements that can cause damage and affect the performance.

A). Rain, Sleet and Snow

If exposed to water, your **MP3R** is susceptible to the damage to electronic or mechanical components. Water can cause electronic malfunction or promote premature corrosion of electrical components and frame.

B). Temperature

Some of the parts of **MP3R** are susceptible to change in temperature. The controller can only operate in temperatures that range between 18 F (-8 C) and 122 F (50 C).

At extreme low temperatures, the batteries may freeze, and your **MERITS MP3R** may not be able to operate. In extreme high temperatures, it may operate at slower speeds due to the controller's safety feature to prevent damage to the motors and other electrical components.



SPECIFICATIONS



MP-3R

MP3R..

Specifications	English
Weight	
Base	82 lbs
	33 lbs
Batteries (NF 22), each	31 lbs
Wheel Size	
Drive Wheels	14 in
Casters	8 in
Anti-tippers	5 in
Seat Dimensions	
Width (adjustable)	18 - 24 in
Backrest Height (adjustable)	19 - 24 in
Depth	18 in
Footrest	
Height adjustable	3 - 7 in
Angle adjustable	0° - 8°
Controller	Integral, Dynamic DL , fully programmable
Charger	On-board or portable, 110 - 220 V, 50 - 60 Hz, 5A
Colors	Cherry, Dark Green
Battery Type (recommended)	2 (two) 22 NF, sealed lead acid or gel cell

ASSEMBLY and DISASSEMBLY

Your MP3R is shipped partially disassembled to maximize its protection during shipping. Please follow the steps described below to assemble it for immediate use:

Tools required for adjustment:

Allen Keys (5, 17 m/m), Wrenches (10, 17 m/m), Screw Drive

Components:

1. Frame base
2. Shroud
3. Seat assembly w/ Seat belt
4. Armrests
5. Integral controller/bracket

Assembly:

Only four components are needed for assembly in the following order.

A). Install Battery & Ignition key connection:

1. Place the batteries to your MP3R. Be sure to place the Positive (+) contact of each battery adjacent to the Negative (-) of the other battery, as shown in (FIG. 1).

 **Note:** The batteries weigh 31 lbs. / each. Please seek help if you cannot lift them by yourself.

2. Connect Positive (+) battery terminals, red, one at a time using screws, washers and nuts provided.
3. Connect Negative (-) battery terminals, black, one at a time.
4. Plug the color-coordinated connectors (batteries and motors).
5. Cover battery post with caps.
6. Connect the ignition cable to the ignition key lock inside the shroud.

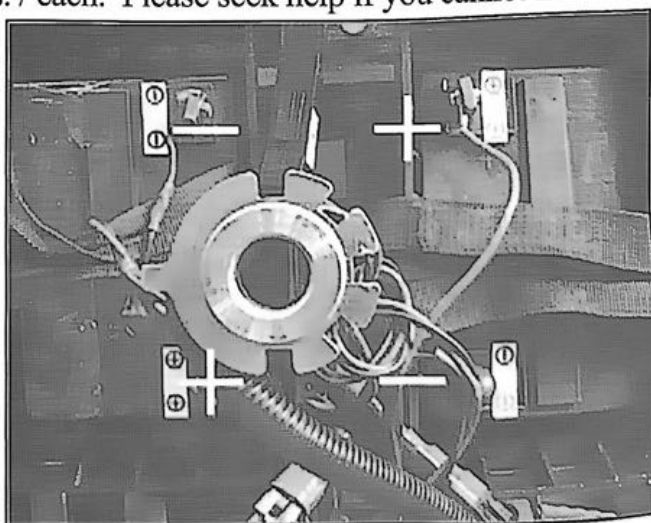


FIG. 1

7. Slide the shroud through the seat pedestal and main (controller) cable.
8. Shroud should also slide through motor free-wheeling levers which should be positioned in engaged mode.
9. Secure shroud to frame.

B.) Installing the Seat:

1. Lift the seat and slide the seat post (bottom of the seat) into the seat pedestal.



Notes: Folding the seat can make it easier to install. (FIG. 2)

2. Be sure to swivel the seat and push down on the seat until it locks in place.

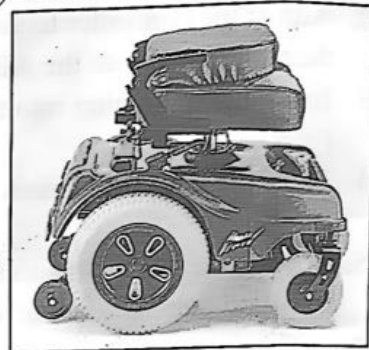


FIG. 2

C.) Inserting the Height and Width adjustable Armrests:

Width

1. Loosen the knobs on the armrest receiver. (FIG. 3)
2. Slide armrest into the horizontal receiver brackets.
3. Select desired width and tighten the knobs.



FIG. 3

Height

1. Locate and loosen the knobs on the vertical armrest holder. (*FIG. 3*)
2. Insert the armrest into the holder.
3. Select desired height and tighten the knobs as hard as possible.

D). Installing the controller:

1. Insert controller bracket tube into the receiver.
(*FIG. 4*)
2. Adjust the controller to your desired length,
then tighten it with the Allen screw.
3. Insert the main plug into the controller socket.
(*FIG. 5*)
4. Place the connected main socket into the Socket-Holster, located under the seat.
(*FIG. 6*)
5. Secure the main cable to the seat using tie-wraps to prevent it from getting in the
way and become damaged. (*FIG. 7*)

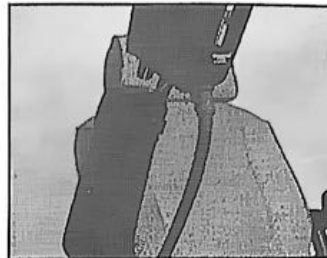


FIG. 4

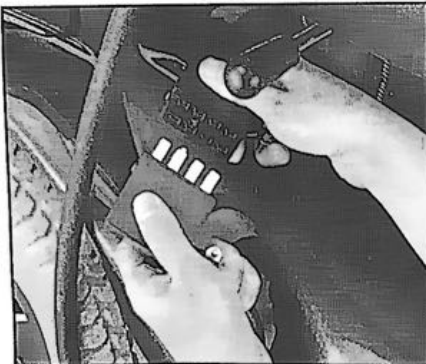


FIG. 5

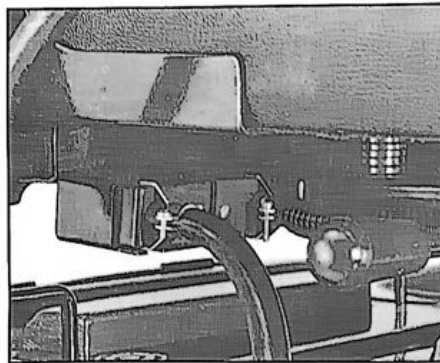


FIG. 6

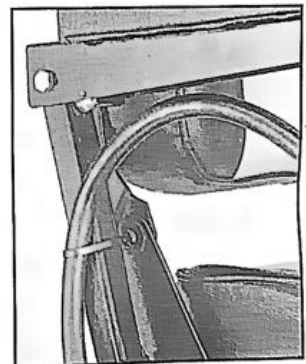


FIG. 7

Notes: Do not tighten the tie-wraps too tight to allow sliding of it when swiveling the seat.

Disassembly:

The disassembly is the reversal of assembly.

A). Seat Removal

1. Turn power off.
2. Make sure that the motors are engaged (not in free-wheel mode).
3. Disconnect main cable from controller.
4. Push the swivel lever forward to unlock it and rotate the seat a little bit to the position where you can feel the seat rotate freely.
5. Pull seat upward.

B). Shroud Removal

After removing the seat

1. Lift shroud straight up.
2. Unplug ignition wires.

C). Battery Removal

i) For transportation:

After removing the seat and shroud,

1. Unplug battery connectors (2, black).
2. Lift batteries one at a time using battery straps.

ii) For replacement:

After removing the seat and shroud,

1. Unplug battery connectors (2, black) to prevent damage to controller.
2. Remove battery post caps.
3. Unscrew positive terminals (+) one at a time.
4. Unscrew negative terminals (-) one at a time.
5. Place battery connectors (wire harnesses), screws, washers and nuts in a safe place.
6. Remove batteries one at a time using battery straps.

ADJUSTMENTS FOR SEATING COMFORT

To maximize your seating comfort, your MP3R lets you adjust the seat width, height and backrest angle; armrest height and tilt angle; and the footrest height and angle. In addition, it also allows you to select your desired controller bracket length.

A). Seat Height

To adjust the seat height .

1. Turn the power off.
2. Disconnect the main cable and remove the seat. (Remove the shroud, if it is necessary).
3. Using a 17 mm wrench, remove the bolt securing the inner seat pedestal. (FIG. 8.1)

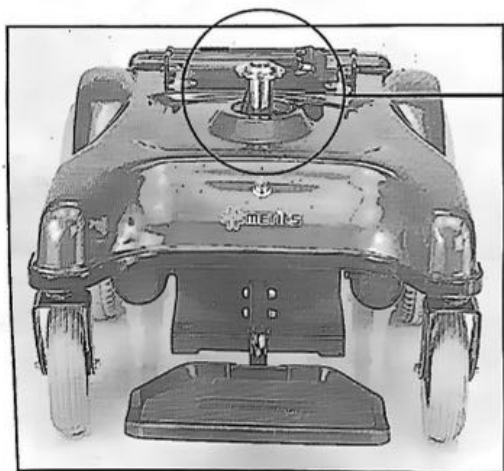


FIG. 8

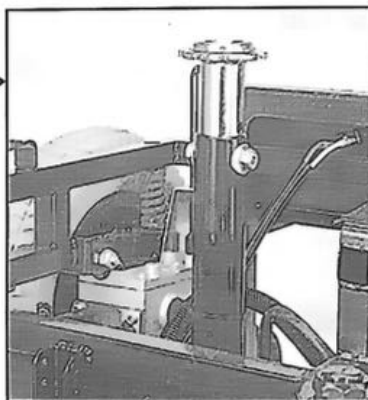


FIG. 8.1

4. Slide the inner pedestal to the desired height and reinsert the bolt.
5. Tighten the bolt and place seat back.

B). Armrest Width, Height, and Angle Adjustment

To adjust the Width

1. Locate the knobs on the armrest holder and turn them counter-clockwise to loosen them. (FIG. 9)
2. Slide the armrest up or down to the desired position.
3. Be sure to retighten them by turning the knobs clockwise as hard as possible.

To adjust the Height

1. Locate the knobs on the armrest holder and turn them counter-clockwise to loosen them. (FIG. 9)
2. Slide the armrest up or down to the desired position.
3. Be sure to retighten them by turning the knobs clockwise as hard as possible.

To adjust the tilt angle

1. Flip up the armrest for easy access.
2. Loosen the hex nut under the armrest pad by turning it counterclockwise with a 5 mm hex wrench. (FIG. 10)
3. Turn the screw counter-clockwise to raise the armrest and clockwise to lower the front of armrest.
4. Be sure to retighten the nut.

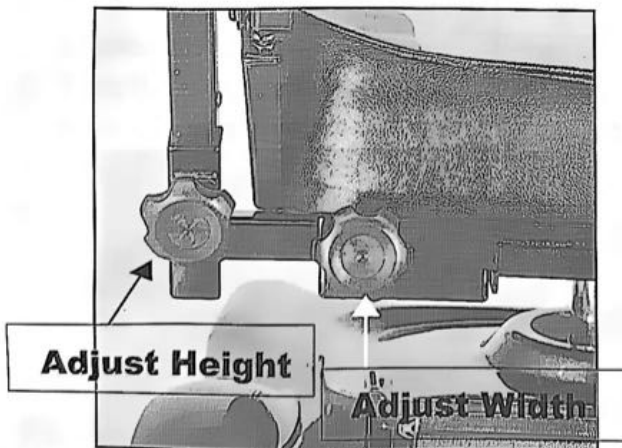


FIG. 9



FIG. 10

C). Foot Platform (*up to 2", 1/2" variance*)

To Adjust Height

After removing the seat and the shroud,

1. Using a 10 mm hex wrench, remove the bolts and nuts. (FIG. 11)
2. Slide the platform to your desired height.
3. Replace the bolts and nuts and be sure to tighten them. (FIG. 11)

To adjust Angle

1. Flip up the footplate for easy access and loosen the nut.
2. With an Allen key, simply turn the bolt counter-clockwise to increase the angle or clockwise to decrease it.
3. Be sure to retighten the nuts. (FIG. 12)

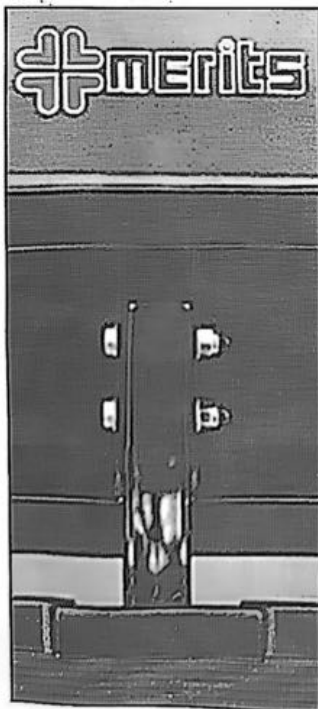


FIG. 11

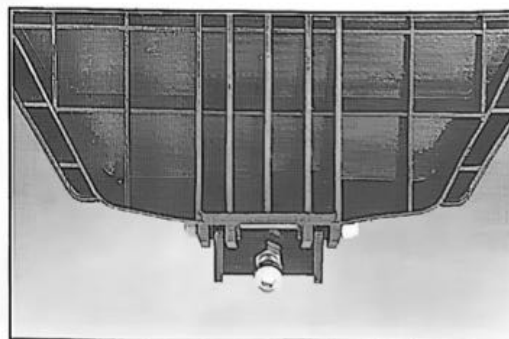
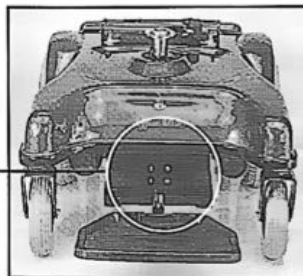


FIG. 12

D). Controller

To adjust the controller length forward or backward

1. Flip up the armrest for easy access.
2. Loosen the bracket bolt with an Allen key.
Slide the controller bracket in or out to your desired length. (FIG. 13)
3. Retighten the bolt.

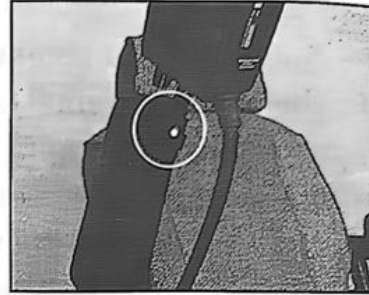


FIG. 13

To reverse controller

1. Disconnect main cord.
2. Remove both sets of armrests, while the controller still is secured on one of armrest.
3. Exchange both armrests.
4. Be sure to tighten the knobs.

E). Backrest Angle

1. Fold down the backrest for easy access.
2. Loosen the nut.
3. Turn the screw clockwise to increase the backrest tilt angle and counter-clockwise to decrease it.
(Fig. 14)
4. Be sure to retighten the nuts.

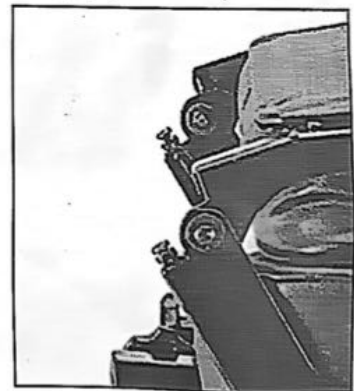


FIG. 14

F). Anti-tipper is adjustable

1. Loosen the bolt by using 6 mm wrench.

2. Adjust the height to one of third position where you desire. (Fig. 15)
3. Be sure to tighten the bolt and nut.



Note: It is not recommended to change the original setting to prevent from tipping over when raising the anti-tipper and from anti-tipper raising the driving wheel on any begin of ramp while lowing the anti-tipper.

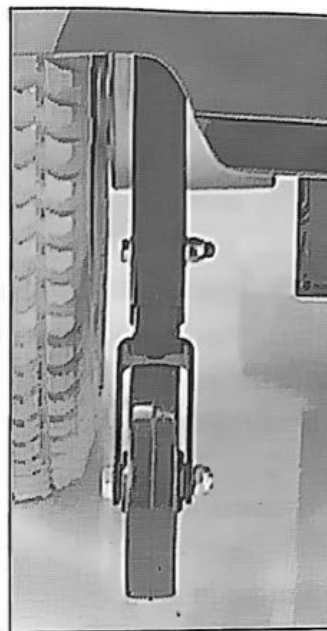


FIG. 15

OPERATING

The MP3R is simple to operate. However, we recommend that you read carefully the following instructions to get familiarized with your new vehicle.

A Word of Caution:

In addition to the power ON/OFF switch that comes standard with the controller, we have included an ignition switch that will allow you to completely disable the operation of the vehicle by removing the ignition key.

Before you turn the power on, always be aware of the environment that surrounds you to select your desired speed. For indoor environments, we recommend that you select the lowest speed setting. For outdoor operation of this vehicle, we recommend that you select a speed that is comfortable for you to control it safely.

The following are the steps and the components required to safely operate your vehicle:

A. Driving:

1. Ignition Switch:

Insert the key in the ignition and turn it clockwise (ON).

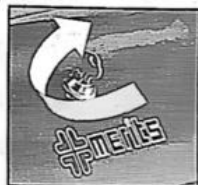


FIG. 16

2. Controller On/Off Switch

Depress the ON/OFF button (I/O) switch, located in front of the joystick to activate your MP3R. The battery gauge lamp will light up to indicate the current charge of your battery. Depressing the ON/OFF button again will deactivate the controller.

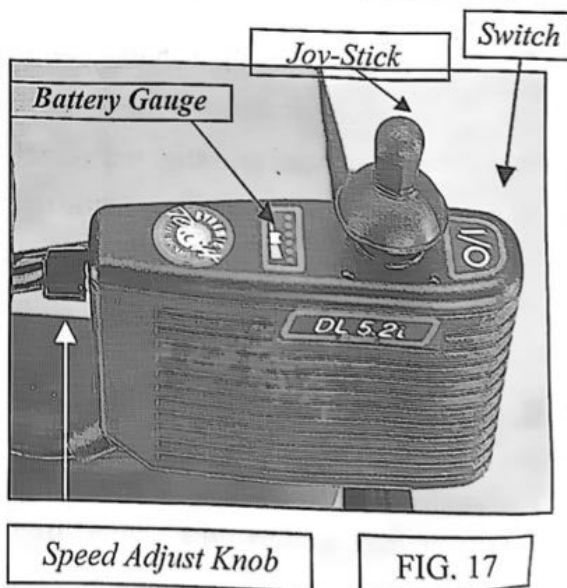


FIG. 17



Notes: If the key is not in the "ON" position when you try to turn the controller ON, the battery charge display will begin flashing, and let you know that the ignition key must be turned on before turning the controller on.

3. Speed Control

The speed control knob allows you to set the forward speed to your desired setting. Turning the knob fully counter-clockwise will set the speed at the slowest setting, while turning the knob fully clockwise will set the speed at the fastest setting. The controller sets the reverse speed, acceleration and deceleration proportionally and automatically for your safety.

4. Joystick

The joystick controls the direction and speed of your vehicle. Pointing the joystick away from the neutral position, center, will move the vehicle in the direction where the joystick is pointing. The farther away forward/backward the joystick is from the neutral position, the faster the vehicle will go. The farther away to the right/left the joystick is pointing, the sharper the turn of the vehicle will be.

To operate the vehicle, gently push the joystick in the direction you want to go. Returning the joystick to its neutral position, center, will reduce the speed and stop the vehicle by automatically applying the brakes (electromagnetic). Gentle operation of the joystick will result in smoother transitions in speed and direction, while sharp operation of the joystick will result in drastic transitions in direction and velocity.

B. Controller Display:

The controller display is a multifunction visual display. It can provide three types of information: a). ON/OFF status, b). battery charge and c). fault diagnostics.

i) ON/OFF Status

When the power is on, the controller's LED will be lit up. If the LED is not lit, the controller is not ON.

ii) Battery Charge

The controller LED is composed by 6 segments, two each red, orange, and green. At full battery charge, all 6 segments of the controller LED are lit continuously. With progressive discharge of batteries, successive segments will extinguish in decanting order. A single lit red LED represents the lowest state-of-charge.

When the battery drops to the low state-of-charge, a warning is given by slowly flashing the left-most read LED on and off. Depending on age and condition of the batteries, you are, then, left with a limited driving range before you must recharge the batteries.

Recommendation:

To ensure a dependable battery charge, we recommend to charge the batteries overnight. It will not only spare you unpleasant situations on route, but will prolong the batteries' service life as well.

iii) Diagnostics

The controller's LED can also let you know about problems with any of the vehicle components and the possible location of them. The following chart lists possible problems indicated by the flashing meter.


Diagnosis	Flash Code Sequence
1. Low Battery Voltage Fault	0 _ 0 _ 0 _
2. High Battery Voltage Fault	00 _ 00 _ 00 _
3. Left Motor (or conntection) fault.	000 _ 000 _
4. Right motor (or connection) fault	0000 _ 0000 _
5. Left or Right Park Brake (or connection) fault	00000 _ 00000 _
6. Controller fault	000000 _ 000000 _
7. Motor stalled or joystick out of neutral time out	0000000 _ 0000000 _

"o" means all 6 LED's flash **on** for a period of 0.2 seconds and off for 0.8 seconds.

"—" means all 6 LED's **off** for a period of two seconds before flashing sequence repeats.

C. Freewheeling:

Because the motors are designed to engage the electromagnetic brakes when the vehicle is not in use or when the power is OFF, they also have a manual feature that allows them to "free-wheel." Free-wheeling is accomplished by turning the free-wheeling levers to the free-wheeling position as shown.(FIG. 17)

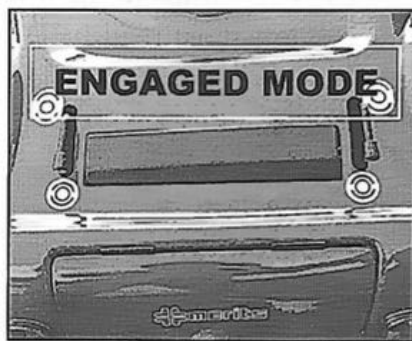

WARNING!

1. Never free-wheel your MP3R on a slope;
2. Never free-wheel the motors while operating your vehicle;
3. Always remember to engage the motors before turning the power back ON.

D. Thermal Rollback:

Your MP3R is equipped with a safe device, thermal rollback circuit. This circuit can monitor the temperature of the controller and motors. In the case of excessive heat of the controller and motors, the controller will decrease the speed of your MP3R to 1/2 of normality on its accord and allow the components to cool down.

Your MP3R will resume its normal speed when the temperature returns to a safe condition.

E. Main Circuit Breaker:

Underneath the rear lid is the charger tray. There is a red button in the hold. (FIG. 19) The main circuit breaker monitors the current draw of battery. The main circuit breaker is a safety feature built in your MP3R for your extra safety. When the batteries and /or the motors are heavily strained, the main circuit breaker might trip for extra protection of your safety. If it happens, please wait for approximately 1 minute. Then, the unit will resume its normal operation if the temperature return to a safe level.

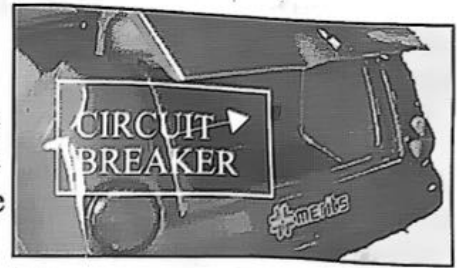


FIG. 19

Usually, the Thermal Rollback Circuit in controller is more sensitive than the Main Circuit Breaker. It is recommended to turn off the power and wait for 5 minutes to allow the components to cool down if you find that you have lost speed or power suddenly.

F). Seat Rotation

Your seat is subject to 90° swivel to each side of your MP3R. (FIG. 20 & 21) The seat lock lever will lock the seat rotation in one of five positions. To rotate the direction of your seat, simply pushing the lock lever forward to unlock and then release it while rotating the seat until it goes to lock in place.



FIG. 20



FIG. 21

BATTERY and CHARGER

A). General Battery Information

This lecture is to introduce the general battery knowledge and information for you. Hopefully, the lecture can provide the useful information for you to allow the battery reaching its best potential performance and life cycle service. In general, the battery should last for years if you take care the battery properly. Otherwise, the battery may die very soon.

Type of Battery Suggested Use

Type: Deep Cycle (Sealed Lead Acid or Gel Cell)
Size: Group 22
Voltage: 12V each
Amperage: 44 (amp hour)

This type of battery is well sealed to prevent leakage and designed without gas when charging and without maintenance. They are recharged by an off-board charging system, on-board charging system, or portable charging system that is for your ease and convenience. For the most efficient, we prefer simultaneous charging of both batteries by the supplied MP3R charger.



Notes: If you insisted to use a different charger, we recommend the battery charger in the 4 to 7 amp ranges. Exceeding the 7 amp will result in damaging heat.

B). Character of this type Battery

Over Discharge

This type of charge is very sensitive to Over-Discharge. And over-discharge can result in failure to recover normal capacity, reduce capacity, or shorten service life.

Retention and Shortage

When the charged battery is left standing for an extended period of time, its

capacity gradually decreases and enters to discharged state. The battery consumes the stored electric energy without releasing it effectively to the circuits. This is called self-discharge. If the shortage is too long, the battery, sometime, cannot recover back its previous capacity.

How does the charge work ?

The charge current is initially large, and decreases towards the end of charging. Therefore, when the battery voltage is low, the charger works very hard to bring up the battery charge. When the battery is fully charged, the charger maintains but does not overcharge the battery.

When do I need to charge the batteries?

It depend on how your usage of your **MP3R**. Also, many factors such as smooth or rough terrain, flat areas or inclines, speed, and other else can affect the battery's consumption. If you used your **MP3R** daily, charging it as soon as you finished using it. The charger will not overcharge your batteries. If you drive your **MP3R** infrequently, you just need to charge it per week for 12 to 14 hours.



Notes: It is recommended to charge the battery per day for over right for the best potential performance and cycle life of your battery.

Store Battery

1. If you do not use your **MP3R** regularly, we recommend maintaining battery animation by charging batteries at least once per week.
2. If you do not plan using your **MP3R** for an extended period time, fully charging the batteries prior to storage and disconnect all connectors on the battery.
3. Store the batteries in a warm environment and avoid temperature extremes.

C). Portable Charging Instructions

Our unique portable charger system is designed for your safety and convenience. This system can provide the easy-to-charge as on-board charge system does while preventing you from electric shock that on-board charge system may confront.

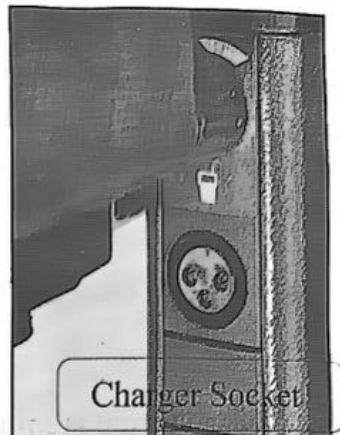


FIG. 22

1. Position your **MP3R** to a standard wall outlet.
2. Turn main power off.
3. Take out the charger from Gumini's charger tray. (Fig. 23)
4. Plug the connector outlet of the charger into the charging socket on the rear of controller.
5. Plug the charger power cord into a standard wall outlet. (FIG. 22)
6. Disconnect the charger power cord from the wall outlet when the battery is full.



FIG. 23

ACCESSORY

There are many accessories provided to your MP3R to make your life more pleasure.


❖ Power Seat

The power seat can be pre-installed by Merits Health Products. 6* infinitive increase can make your life even more independent.




❖ Swing away Legrest

Your MP3R can be fitted to a swing away legrest very easily to provide your feet more comfortable resting.

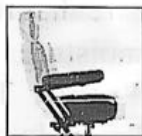


❖ Slide Seat Base

The Slide Seat Base can make your seating position more comfortable. You can select one of eight positions of forward or backward for your unique needs.

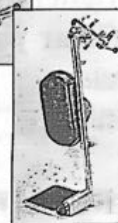


❖ Crutch Holder and Oxygen Holder

Every friendly design Holders, no tolls required to install the holders.

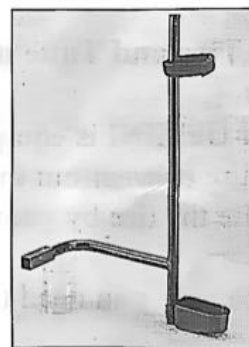

❖ Solid Tires

The wheels of MP3R are standard pneumatic (air filled) but can be ordered as solid. The solid tires are the same size and tread but filled with a foam insert. Solid tires require no maintenance and eliminate the possibility of a flat tire.

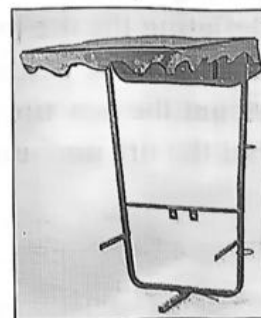

❖ CANOPY

It is designed to protect you from the sun and rain.

MAINTENANCE and REPAIR

Your **MP3R** is well designed for perfection. You just need a very simple maintenance.

A). MAINTANCE:

- **Tire Pressure:** Be sure to maintain pressure of tire between 30-37 psi.
- **Tire Tread:** If the tire tread was less than $1/32^*$, please exchange your tire as soon as possible.
- **Motor Brushes:** Please allow your dealer to inspect the motor brushes for every six months.
- **Joy-Stick:** Be certain to keep the controller from moisture.
- **Battery Terminal Connections:** Be sure that they are not corroded. Battery is needed to keep away from moisture.
- **Clear your **MP3R** with a Damp Wrap:** The ABS shroud has a clear coat that is very easy to clear. Please do not use water to wash the **MP3R**.

B). REPAIR:

◇ Tire and Tube of drive wheel:

Your **MP3R** is equipped with a 14" drive wheel. If you encounter with the flat tire or wear out the tire, you may need to exchange the tire or tube or inflate the tire by yourself.

In the case you need to replace the new tire or tube for the driving wheel:

1. Remove the driving wheel with hexagonal socket wrench, 17 m/m by turning counter-clockwise.
2. Deflating the tire pressure as much as possible.
3. Dismount the wearied tire and tube from the rim.
4. Mount the new tire and / or tube on the rim.
5. Flat the tire pressure to your desire pressure (30-37 psi).

✧ **Tire and Tube of Caster:**

1. Remove the wheel with the wrench with the hexagonal wrench.
2. Deflating the tire pressure as much as possible.
3. Dismount the wearied tire and / or tube from the rim.
4. Mount the new tire and / or tube to the rim.



Notes: If you encounter any technical problems, it is recommended to check with your local dealer.

✧ **Inflate Tire (Emergency) :**

1. Take out the adapter tube from the pump and connect to the pump.
2. Insert the adapter to the tire valve.
3. Start pumping.

LIMITED WARRANTY

Life Time Warranty:

We will repair or replace for the original purchase, free of charge, of following part found upon examination by any authorized service center on by us to be defected in material and./ or workmanship.

- Main Frame

One Year Limited Warranty:

For one year from purchase **MP3R** we will repair or replace for the original purchase, free of charge, any of following parts found upon examination by any authorized service center on by us to be defected in material and / or workmanship.

- Electronic Controller
- Drive Train Components

The one year warranty excludes the

- Motor brushes: The motor brushes are wear items.
- Brake Pad: Brake pad is a war item.

Six Month Limited Warranty:

For six months from date of purchase **MP3R** we will replace for the original purchase, free of charge, any of following parts found upon examination by any authorized service or by us to be defected in materials and / or workmanship.

- Plastic parts
- Rubber parts
- Bearings
- Other parts not specifically identified above.

Warrants Exclusions:

- ABS plastic shrouds

- Batteries (battery are warranted by the battery manufacture)
- Tires and Tubes
- Damage caused by:
 1. Improper installation on the parts of the individual owner
 2. Abuse, Misuse, Negligence, or Accident.
 3. Commercial use
 4. Circumstance beyond the control of us

We sole liability extends to repairs or replacement of parts or components as previously specified. The foregoing is in lieu of all warranties, expressed, implied or statutory, whether written or oral. We should not be liable for any consequential or incidental damages whatever.